# Drive DeVilbiss Moving & Handling











www.drivedevilbiss.co.uk



## Contents

A wide range of equipment enabling you to turn, lift and transfer your clients in a safe environment.

Able Assist Patient Transfer Aids	04
Lifting/Turning Aids	06
Moving & Handling Accessories	10
Slide Sheets	11
Serena Pillow Lift	12
WendyLett Turning System	13
Easy Slide Sheets	14
Bed Management System	19
Casa Lift Hoist	21
Samsoft Hoists	22
Parklite Aluminium Hoists	24
Novaltis Lift Aid	26
Hermes Hoist	28
Portable Bed Hoist	29
Slings	31
Riser Recliner Armchairs	34
Integra Shell Chair	36
Mountway Solo Toilet Lift	38
AquaJoy Premier Plus	40
Mountway Neptune	41
Bathroom Accessories	42

#### Able Assist Patient Transfer Aids

The Able Assist Patient Transfer Aid provides a safe and easy solution for a Carer wishing to transfer a client. A standard version is available (PT003), an adjustable leg version (PT004) and a folding version (PT005).

The Able Assist enables Users to participate in transfers, improving muscle tone. It is ideal for transferring Users from room to room, increasing the User's independence.

For those Users who can raise from a seated position to a standing position, the Able Assist is ideal. The Carer will rotate the seat parts, the User stands on the footplate provided. The Carer will then lower the seat flaps and the User can be seated, having their legs supported by the knee brace, and hold on to the horizontal bar in front of them, making them feel secure. Model PT004 which is easy to operate by pedal, adjustable legs is ideal for getting close up to a Riser Recliner Armchair.

The transfer aid has 4 wheels, 2 of which have brakes allowing this aid to be manoeuvred easily from room to room. The narrow width of this Aid enables the Carer to manoeuvre easily through standard doorways. The Able Assist can be wheeled directly up to a toilet, with wheels positioned either side of the W/C, offering the User more dignity. To sit down, the above process is reversed.

DESCRIPTION	CODE	
Able Assist Transfer Aid - Standard	PT003	
Able Assist Transfer Aid - Leg Adjustable	PT004	
Able Assist Transfer Aid - Folding	PT005	
Transfer Sling	4200BE	

#### PT003 Specifications

- » Width: 61cm (24")
- » Height: 97cm (38")
- » Length: 86cm (34")
- » Knee Pad Width: 39cm (15")
- » Knee Pad Height: 26cm (10") » Knee Pad Height from Base:
- 17 41cm (6.5 16")
- Handle Height: 85cm (33.5")
- » Seat Height to Bottom: 62cm (24.5")

at	Height	to	ton:	67cm	(26")	

- Seat Height to top: 67cm (26")
   Chassis Length: 75cm (29.5")
- Chassis Width (Front): 46cm (18")
- Chassis Width (Back): 61cm (24")
- Chassis Width (Internal): 35cm (14")
- Brake Castors (Back): 7.5cm (3")
- Max User Weight: 185kg (29st)

#### **PT004** Specifications

- » Width: 53cm (21")
- Height: 106cm (42")
- Length: 82cm (32")
- » Knee Pad Width: 39cm (15")
- » Knee Pad Height: 26cm (10")
- » Knee Pad Height from Base:
- 25 49cm (10 19") Handle Height: 90cm (35.5")
- » Seat Height to Bottom:
- 62cm (24.5")
- » Seat Height to top: 67cm (26")
- Chassis Length: 87cm (34")

- Chassis Width (Front): 52cm (20.5") > Chassis Width (Back):
- 50cm (19.5") » Chassis Width (Internal):
- 49 76cm (19 30") > Plain Castors (Front):
- 7.5cm (3") » Brake Castors (Back):
- 7.5cm (3") » Max User Weight:
- 185kg (29st)

- » Plain Castors (Front): 7.5cm (3")

- Width Open: 61cm (24")
- » Width Folded: 37cm (14.5")
- » Height: 97.5cm (38")
- » Length: 78cm (31")
- » Knee Pad Width: 36cm (14")
- » Knee Pad Height: 25cm (10")
- » Knee Pad Height from Base: 19 - 43cm (7.5 - 17")
- » Handle Height: 87cm (34")
- » Seat Height to Bottom:
- 59cm (23")

- Seat Height to top: 77cm (30")
- Chassis Length: 79cm (31")
- » Chassis Width (Front): 46cm (18")
- » Chassis Width (Back): 61cm (24")
- » Chassis Width (Internal): 34.5cm (13.5")
- Plain Castors (Front): 7.5cm (3")
- » Brake Castors (Back): 7.5cm (3")
- » Max User Weight: 185kg (29st)

## Lifting / Turning Aids

#### Orbital Transfer Platform

The Orbital Transfer Platform assists a person who is weight bearing and has good upper body strength to transfer from one seated position to another.

- Specifications
- Width: 66cm (26")
  Height to Front Handles: 76cm (30")
- Height to Top Bar: 99cm (39")
- Depth: 69cm (27")
- Max User Weight:
- 222kg (35st)

CODE	DESCRIPTION
800	Orbital Transfer Platform without Knee Support
805	Orbital Transfer Platform with Knee Support
800FA	Orbital transfer platform with Forearm Platform
806PUSUPPORT	Knee Support Only



#### Ablestand Patient Turner

The Ablestand Patient Turner helps a Carer or attendant to transfer a User to another seat, bed, wheelchair, toilet or commode by standing and turning. Easily dismantled/ reassembled without tools for transportation and storage.

#### Specifications

- » Width: 49cm (19")
- Height: 80 126cm (31.5 50")
   Depth: 45cm (18")
- Maximum User Weight:
- 150kg (24st)

CODE	DESCRIPTION
PT002	Ablestand Patient Turner



#### Ambiturn Sit to Stand Transfer Aid

The Ambiturn Sit to Stand Transfer Aid assists with pivoting transfers between different seat surfaces and can transport Users short distances.

#### Specification

- » Width: 63cm (25")
- » Height: 118cm (46.5")
- Depth: 57.5cm (22.5")
  Max User Weight:
- Max User Wei
   150kg (24st)

CODE	DESCRIPTION
PT001	Ambiturn Sit to Stand Transfer Aid

#### **Turner Pro**

The Turner PRO is a turn aid with a functional design that offers safe patient turning with standing support. The only adjustment necessary is the height of the leg support and the handle. Thereafter, the transfer can be performed safely whilst the caregiver and user can maintain eye contact and communication. The orange handle provides a strong, high contrasting colour which can help users with impaired vision or dementia. The small footplate makes it the perfect choice of transfer aid when transferring in tight spaces.

#### Specification

- » Width: 38cm (15")
- » Height: 78 120cm (30 49")
- » Depth: 45cm (18")
- Maximum User Weight: 2004rg (21 Ect)
- 200kg (31.5st)

CODE	DESCRIPTION
16090106	Turner Pro















## Lifting / Turning Aids

#### Ambiturn Sling

For additional safety and peace of mind a transfer sling can be used with the Ambiturn. The transfer sling has a Velcro strap, in addition to a buckle which locks in place, the strapping of which fits over the handle bar of the Ambiturn transfer aid, enabling the User to feel more secure.

CODE	DESCRIPTION
PTBXXS	Extra Extra Small
PTBXS	Extra Small
PTBS	Small
PTBM	Medium
PTBL	Large
PTBXL	Extra Large



#### Ambiturn and Ablestand Deluxe Feeder Sling

The Deluxe one Carer Feeder Ambiturn/Ablestand Sling enables the Carer to assist the User to a standing position, using the Ambiturn or Able Stand. Carer assists the user to stand by feeding the strap through the belt. Available in three sizes - Small, Medium and Large

CODE	DESCRIPTION
4401BE	Small
4402BE	Medium
4403BE	Large



#### Turn Table

This strong and lightweight turntable provides easy rotation. Swivels 360° for smooth and easy moving in any direction. It enables the User to be transferred more safely from chair to chair, from chair to bed etc.

- » Diameter: 40cm (16")
- Weight: 1kg (2lb)
   Maximum User Weight: 115kg (18st)

CODE	DESCRIPTION
HCA013	Turn Table

#### Swivel Seat

Developed to aid people who have mobility problems to turn on a number of different types of seats. Top disc rotates a full 360° over bottom disc, allowing smooth two way turning action. Suitable for many applications including dining chairs or car seats. Slip resistant underside enables secure positioning on seat. Fleece cover is warm and comfortable and removable for machine washing.

- » Diameter: 40cm (16")
- Maximum User Weight: 115kg (18st)

CODE	DESCRIPTION
SWIV001RT	Swivel Seat





#### Moving & Handling Accessories



**Easy Stool** is a versatile piece of equipment which is used to aid the repositioning of a User in a chair, in to and out of bed or on a toilet. Other applications include facilitating sling fitting, wound dressing, tying shoe laces and raising a User's knees while seated on a toilet.

**Oval Slide Boards** is a sliding board which is used for supine transfers from one bed or stretcher to another but also for positioning higher up in the bed. It can also be used to assist placement of a lifting sling both in the recumbent and in seating positions.



**Easy Gloves** are mits which can be used in a multitude of moving and handling procedures to reduce friction and shear. The easy gloves can be used to release skin loading when profiling in a bed, applying slings or inserting easy slide boards.

#### Leg Lifter is a

comfortable padded leg lifter or handling sling that can be used for transfers or the re-positioning of patient's legs safely on or off a leg rest.

#### **Slide Sheets**

Slide Sheets allow basic handling without the need to lift the patient. They are invaluable when moving a patient on or off a bed, for turning the patient in bed and as an aid to sitting the patient up in bed.

Tubular slide sheets allows 3 basic patient handling tasks without the need to lift and without the risk of shearing to the User's skin.

(1) Turning a patient

(2) Sitting a patient up in bed

(3) Moving a patient on/off bed

The Slide Sheets are made from strong nylon material which allows friction free movements on both sides of the Sheet. Slide Sheets reduce manual handling effort and strain, minimizing the risk of back and other injuries to Carers. They reduce skin tears and bruising and promote patient comfort and dignity. Each Sheet is provided with a hanging loop for storage or handle loops, which are colour coded for size.

CODE	DESCRIPTION	COLOUR LOOP
4301BE	Tubular Slide Mini Reusable (60cm x 41cm (23.5" x 16"))	Blue
4302BE	Tubular Slide Standard Reusable (122cm x 71cm (48" x 28"))	Black
4303BE	Tubular Slide Sheet Long Reusable (145cm x 71cm (57" x 28"))	Green
4304BE	Tubular Slide Sheet Extra Long Reusable (200cm x 71cm (79" x 28"))	Yellow
4305BE	Tubular Slide Sheet Maxi Reusable (200cm x 140cm (79" x 55"))	Red

CODE	DESCRIPTION
5000BE	Easy Stool
5001BE	Easy Gloves (Supplied in Pairs)
5002BE	Oval Slide Boards 46 x 19cm (Supplied in Pairs)
5003BE	Oval Slide Boards 85 x 39cm (Supplied in Pairs)
5004BE	Leg Lifter

Handle loops are colour coded for size

#### Serena Pillow Lift

The Serena Pillow Lift is an electrically powered backrest designed to assist Users with limited upper body strength to alter their position from lying to sitting in bed.

- Also allows Users with breathing or cardiac problems to sleep fully supported in a semi elevated position
- In its raised position, it offers support for reading, eating or watching TV and may also assist with getting in and out of bed
- Fully portable and can be easily fitted either side of any style of bed - slatted, metal framed or divan
- » Simply placed on top of mattress
- Compact frame and ultra quiet operation allows partner to sleep in same bed uninterrupted
- Sturdy, lacquered plywood frame blends discreetly into any bedroom decor and is easily wiped clean with non abrasive household cleaner
- Central actuator is unobtrusive and provides a smooth, stable lifting action
- » Soft foam mattress made from combustion modified polyether foam offers comfort for long periods in bed
- » Dartex mattress cover is waterproof and fully machine washable for hygiene
- Soft touch hand control can be effortlessly operated by arthritic hands
- » Maximum degree of lift: 80°

Soft touch hand control can be effortlessly operated by arthritic hands

CODE	DESCRIPTION	FRAME WIDTH	FRAME LENGTH	LENGTH WITH MATTRESS	MAX HEIGHT UP	MAX HEIGHT DOWN	WEIGHT WITH MATTRESS	WEIGHT CAPACITY
SERENAS	Serena with Short Mattress	71cm (28″)	72cm (28″)	73cm (29″)	66cm (26")	15cm (6″)	10kg (23lb)	135kg (22st)
SERENAM	Serena with Standard Mattress	71cm (28″)	72cm (28″)	124cm (49″)	66cm (26")	15cm (6″)	11kg (25lb)	135kg (22st)
SERENAL	Serena with Long Mattress	71cm (28″)	72cm (28″)	183cm (72")	66cm (26")	15cm (6″)	12kg (26lb)	135kg (22st)

### WendyLett Turning System

The WendyLett Turning System has been developed to facilitate repositioning in bed. The system consists of high and low friction material working together in the design of the fabric: the weave is smooth in one direction only.

- WendyLett's components can be washed at up to 90°C for efficient infection control
- Can remain underneath the User as bed linen compared to traditional slide sheet systems that must be removed
- » Pattern of the fabric (striped or check) indicates direction of slide
- » Anti-slide material at edges for greater safety
- Choice of three styles/sizes of base sheet (midi, maxi or fitted) depends on risk assessment of the User

NAME	PART NO	SIZE	SLIDING AREA
Base Sheet Mini	R1629	100 x 200cm	60cm (23.5″)
Base Sheet Midi	R1634	140 x 200cm	60cm (23.5″)
Base Sheet Maxi	R1639	200 x 200cm	70cm (27.5″)
Base Sheet Fitted	R1641	200 x 90cm	70cm (27.5")
Base Sheet Fitted	R1651	200 x 105cm	80cm (31.5")
Base Sheet Maxi	R1669	200 x 220cm	100cm (39")
2WAY Mini	R1631	100 x 200cm	n/a
2WAY Midi	R1636	140 x 200cm	n/a
4WAY Midi	R1637	140 x 200cm	n/a
4WAY Maxi	R1647	200 x 200cm	n/a
2WAY Mini Complete System (R1629 and R1631)	R1630	100 x 200cm	60cm (23.5")
2WAY Midi Complete System (R1634 and R1636)	R1635	140 x 200cm	60cm (23.5")
2WAY Maxi Complete System (R1639 and R1636)	R1640	200 x 200cm	70cm (27.5″)
2WAY Midi Fitted Complete System (R1641 and R1636)	R1642	200 x 90cm	70cm (27.5″)
4WAY Midi Complete System (R1634 and R1637)	R1638	140 x 200cm	60cm (23.5")
4WAY Maxi Complete System (R1639 and R1637)	R1643	140 x 200cm	70cm (27.5″)
4WAY Midi Fitted Complete System (R1641 and R1637)	R1644	200 x 90cm	70cm (27.5")
4WAY Maxi Fitted Complete System (R1641 and R1647)	R1645	200 x 90cm	70cm (27.5″)

#### WendyLett Base Sheet

- Assists turning in bed for Users with mild to moderate difficulty in turning independently
- Sliding is lateral in the direction of the stripes
- Slide sheet has a clear pattern of wide grey stripes on a white background
- Excellent for sufferers of MS, Parkinsons, etc.
- Offers many cost benefits for service provider and User



#### WendyLett2WAY System

- Assists the Carer with lateral move and fine positioning of the User when turning unaided is not possible
- 2 way system prevents the User sliding down the bed
- Slide sheet has a clear pattern of wide grey stripes on a white background
- User can be turned laterally in the direction of the stripes when placed on top of the WendyLett base sheet



#### WendyLett4WAY System (Multi-Directional)

- Enables repositioning in multiple directions (up/down/sideways) when placed on top of the WendyLett base sheet
- Allows easy turning, repositioning and recumbent positioning of the User in this combination
- Brake system can be created by tucking the draw sheet under the mattress and the User's buttocks & shoulders ensuring User security
- Offers easy and safe positioning of users who have a one-sided weakness



### Easy Slide Base Sheets (Single)


## \*All dimensions shown are in cm's


200cm


### Easy Slide Base Sheets (Double)


The easy-slide sheet is a two-way sheet that is designed to aid people with impaired mobility to turn independently. The base sheet has a low friction panel which can cover the entire bed length or 3/4 of the bed length for users who are able to push themselves with their feet.

The easy-slide base sheet incorporates high friction side panels to stop the user from sliding too far across the bed. The easy-slide base sheet is also ideal for women who are heavily pregnant as it reduces the friction and shearing on the hips when attempting to turn during the night.

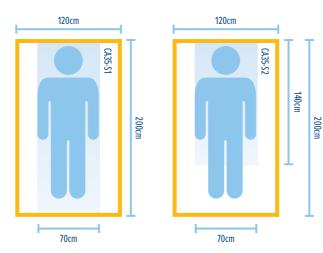

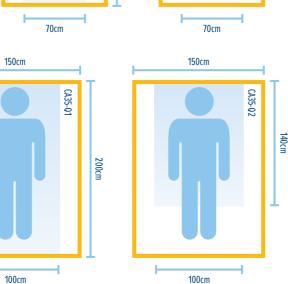

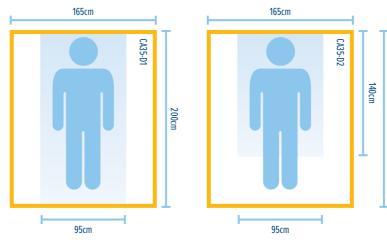

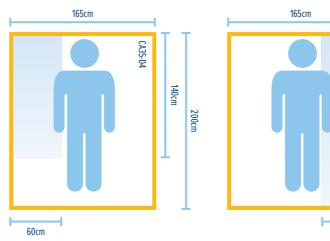


CODE	DESCRIPTION
CA35-S1	Easyslide Base Sheet - 120 x 200cm - Full length satin centre (70 x 200cm)
CA35-S2	Easyslide Base Sheet - 120 x 200cm - 3/4 length satin centre (70 x 140cm)
CA35-Q1	Easyslide Base Sheet - 150 x 200cm - Full length satin centre (100 x 200cm)
CA35-Q2	Easyslide Base Sheet - 150 x 200cm - 3/4 length satin centre (100 x 140cm)

1

CA35-D5

60cm

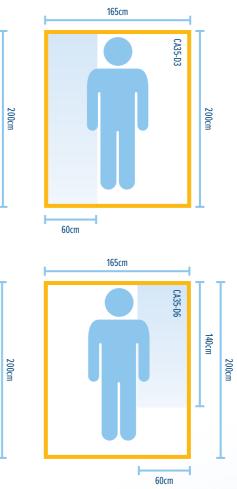

CODE	DESCRIPTION
CA35-D1	Easyslide Base Sheet - 165 x 200cm - Full length satin centre (95 x 200cm)
CA35-D2	Easyslide Base Sheet - 165 x 200cm - 3/4 length satin centre (95 x 140cm)
CA35-D3	Easyslide Base Sheet - 165 x 200cm - Full length satin right (60 x 200cm)
CA35-D4	Easyslide Base Sheet - 165 x 200cm - 3/4 length satin right (60 x 140cm)
CA35-D5	Easyslide Base Sheet - 165 x 200cm - Full length satin left (60 x 200cm)
CA35-D6	Easyslide Base Sheet - 165 x 200cm - 3/4 length satin left (60 x 140cm)

### **Easy Slide Base Sheets (King Size)**

#### \*All dimensions shown are in cm's

180cm

100cm

180cm

CA35-K3

CA35-K6

100cm


140


## **Easy Slide Top Sheets (Single)**

The easy-slide top sheets system designed to work in conjunction with the easy-slide base sheet.


When applied to the easy-slide base sheet, the system allows the user to be repositioned in bed in all directions, due to the under panel of low friction material. The system can be used with either one or two carers.

The system also allows a safer turn of the User as it allows the turn to be compiled in the centre of the bed, which reduces the cases where the User is turned on their side but is too close to the edge of the bed.

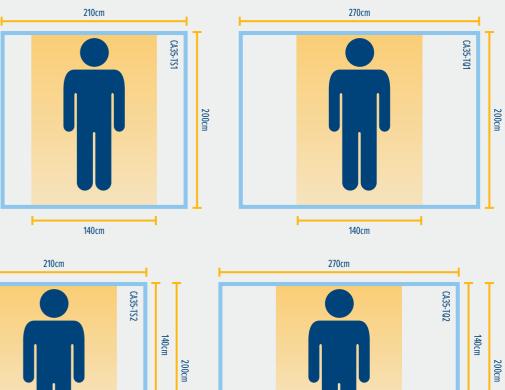

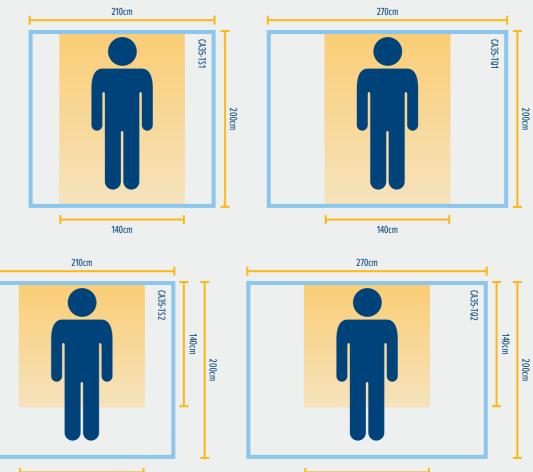

140cm

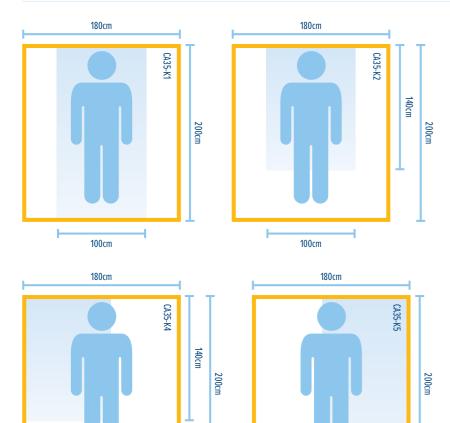



100cm

100cm

CODE	DESCRIPTION
CA35-TS1	Easyslide Top Sheet - 210 x200cm - Full length (140 x 200cm)
CA35-TS2	Easyslide Top Sheet - 210 x 200cm - 3/4 length (140 x 140cm)
CA35-TQ1	Easyslide Top Sheet - 270 x 200cm - Full length (140 x 200cm)
CA35-TQ2	Easyslide Top Sheet - 270 x 200cm - 3/4 length (140 x 140cm)

140cm

### Easy Slide Top Sheets (Double & King Size)

\*All dimensions shown are in cm's


## Bed Management System


The Freeway Bed Management System provides an easy to use, low friction solution for patients who require frequent repositioning.


The Freeway Bed Management System consists of two separate units, the system has been designed to provide comfort for patients and can help to reduce incidences of pressure ulcers. It also provides ease of use by carers and forms an ideal bed solution as part of a single handed care package. Can be used on pressure care mattresses and can also be used with a hoist to assist turning. Dimensions 140 x 150cm. Maximum user weight 200Kg.




CODE	DESCRIPTION
2002AACBMS-COM	Freeway Bed Management System Complete
2002ACCBMS-BSS	Bottom sheet only
2002ACCBMS-QS	Quilted cotton top sheet only

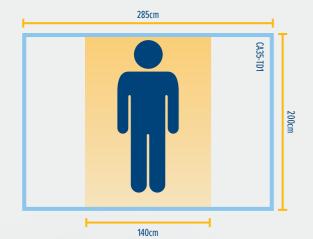

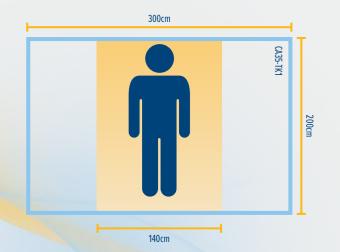

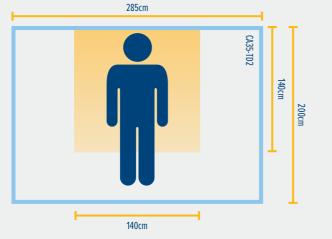

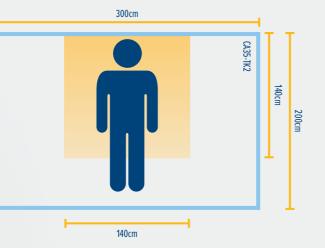


CODE	DESCRIPTION
CA35-TD1	Easyslide Top Sheet - 285 x 200cm - Full length (140cm x 200cm)
CA35-TD2	Easyslide Top Sheet - 285 x200cm - 3/4 length (140 x 140cm)
CA35-TK1	Easyslide Top Sheet - 300 x 200cm - Full length (140 x 200cm)
CA35-TK2	Easyslide Top Sheet - 300 x 200cm 3/4 length (140 x 140cm)



Breathable waterproof quilted cotton top sheet.

Water esistant under sheet with wide non slip sides

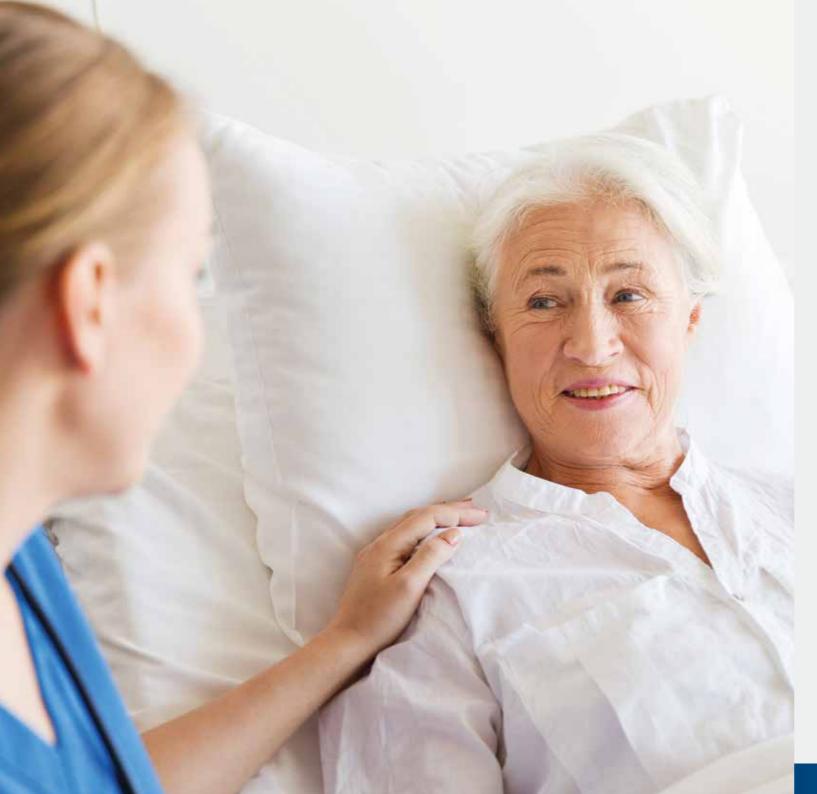

### Casa Lift Hoist

The Casa Lift Hoist 18001BE is a lightweight and versatile lifting hoist that is easy to move into position and allows Users to achieve a lower lifting position even from the floor. This electrically operated hoist has a very high lifting range at a maximum of 170cm (67"), making it one of the highest lifting hoists on the market.

- High lift range of 170cm (67")
   Maximum User Weight Limit 180kg (28st)
   Quality Linak control system with electric leg opener
- » 50,000 lift cycle
- > 50,000 Int cycle
   > Low clearance, easy access underneath beds
   > V shaped spreader bar offers greater head clearance
   > Compact, easy access through doorways
   > Ease of sling attachment
   > Optional counter / service indicator available
   > Optional counter / service and indicator available

- Optional weighing scales available
   3 years warranty frame, 5 years warranty control system, 12 months warranty on battery



DESCRIPTION	CODE
Casa Lift Hoist	18001BE
Casa Lift Hoist complete with Counter / Service Indicator	18001CBE

LIFTING RANGE (BOOM POSITION LOW-HIGH)	CHASSIS LENGTH	CHASSIS WIDTH (SPREAD)	CHASSIS WIDTH (PARALLEL)	CHASSIS HEIGHT (7.5CM)	CHASSIS HEIGHT (10CM)	WEIGHT CAPACITY
44 - 170cm (17 - 67″)	115cm (45″)	112cm (44″)	59cm (23″)	10cm (4″)	13.5cm (5.5″)	180kg (28st)



### Samsoft 150

#### Samsoft 175

The Samsoft Mini's compact design enables it to be manoeuvred easily in small spaces and through doorways. This electrically operated hoist has manual pedals, enabling the chassis legs to be opened with ease and wheeled around a chair.

- Compact design enables it to be manoeuvred easily in small spaces and through doorways
- » Can be vertically stored away when not in use
- Electrically operated hoist has manual pedals, enabling the chassis legs to be opened with ease and wheeled around a chair
- » Quality Dewert control system
- » On / Off board charging system
- » Low clearance, easy access underneath beds
- » Charger holder
- » Colour coded brake castors
- » Comprehensive range of slings available

CODE	DESCRIPTION	LIFTING RANGE (BOOM POSITION LOW-HIGH)	CHASSIS LENGTH	CHASSIS WIDTH (OVERALL)	CHASSIS WIDTH (INTERNAL)	CHASSIS HEIGHT	WEIGHT CAPACITY
18009BE	Samsoft 150 Hoist	64 - 170cm (25 - 67″)	89cm (35″)	56cm (22″)	45 - 79cm (18 - 31″)	9cm (3.5″)	150kg (24st)

CODE	DESCRIPTION	LIFTING RANGE (BOOM POSITION LOW-HIGH)	CHASSIS LENGTH	CHASSIS WIDTH (OVERALL)	CHASSIS WIDTH (INTERNAL)	CHASSIS HEIGHT	WEIGHT CAPACITY
18010BE	Samsoft 175 Hoist	38 - 160cm (15 - 63″)	112cm (44″)	61cm (24″)	50 - 94cm (20 - 37″)	9cm (3.5″)	175kg (27.5st)

Folding hoist can be vertically stored away when not in use

The Samsoft 175 has a greater weight capacity than the Samsoft 150 at 175 Kg (27.5st) and can be vertically stored away when not in use. This electrically operated hoist has manual pedals, enabling the chassis legs to be opened with ease and wheeled around a chair.

- » Higher weight capacity than the Samsoft Mini at 175 Kg
- » Folding hoist can be vertically stored away when not in use
- Electrically operated hoist has manual pedals, enabling the chassis legs to be opened with ease and
- wheeled around a chair
- » Quality Dewert control system
- » On / Off board charging system
- » Low clearance, easy access underneath beds
- » Charger holder
- » Colour coded brake castors
- » Comprehensive range of sling available

#### **Parklite Aluminium Hoists**

180

The Parklite range of hoists are manufactured in lightweight aluminium. Electrically operated user functions allow a single carer to perform safe and secure patient lifts without the risk of injury to their neck or shoulders. The intelligent service diagnostic system with visual LCD display is built in to the control box to make servicing ever easier.

> The cycle counter records every lift cycle when the hoist is loaded, whilst the unique work counter records every second that the actuator is working regardless of whether the lift is loaded or unloaded. Usage information is displayed visually making identifying the age and condition of the actuator very easy. The Parklite range of hoists have been specifically designed for the Acute Sector, satisfying the highest standards for safety and reliability.

320

CODE	DESCRIPTION	LIFTING RANGE (BOOM POSITION LOW-HIGH)	CHASSIS LENGTH	CHASSIS WIDTH (OVERALL)	CHASSIS WIDTH (INTERNAL)	CHASSIS HEIGHT	WEIGHT CAPACITY
4001010001	Parklite 180 Hoist	44 - 173cm (17 - 68″)	114cm (45″)	61cm (24″)	52 - 99cm (20 - 39″)	10.5cm (4")	180kg (28st)
4001010003	Parklite 320 Hoist	53 - 182cm (21 - 71″)	150cm (59″)	79cm (31″)	65 - 105cm (25 - 41″)	15.7cm (6″)	320kg (50st)

#### **Parklite Standing Raiser**

The Parklite stand aid hoist is manufactured from lightweight aluminium. Electrically operated user functions allow a single carer to raise the User from a seated to standing position, without the risk of injury to their neck or shoulders.

The intelligent service diagnostic system with visual LCD display is built in to the control box to make servicing ever easier. The cycle counter records every lift cycle when the hoist is loaded, whilst the unique work counter records every second that the actuator is working regardless of whether the lift is loaded or unloaded. Usage information is displayed visually identifying the age and condition of the actuator.

## Transport

Padded polyester material. Velcro and buckle fastening. Coloured loops. Legs for added security

\*Slings fit all Parklite and Novaltis Hoists

CODE		DESCRIPTION					
3901TBE		3901TBE Small Transport Sling for Stand Aid Hoist					
3902TBE		3902TBE Medium Transport Sling for Stand Aid Hoist					
3903TBE		3903TBE Large Transport Sling for Stand Aid Hoist					
			1				
CODE	DESCRIPTION	LIFTING RANGE (BOOM POSITION LOW-HIGH)	CHASSIS LENGTH	ł	H CHASSIS WIDTH (OVERALL)	H CHASSIS WIDTH (OVERALL) CHASSIS WIDTH (INTERNAL)	H CHASSIS WIDTH (OVERALL) CHASSIS WIDTH (INTERNAL) CHASSIS HEIGHT
4001010002	Parklite 200 Hoist	94 - 144cm (37 - 56″)	95cm (37")		57cm (22″)	57cm (22") 46 - 71cm (18 - 28")	57cm (22") 46 - 71cm (18 - 28") 13.5cm (5")



### **Novaltis Lift Aid**

The Drive Novaltis enables Users to be quickly and efficiently transferred with ease. The Novaltis requires the User to participate with the lift/transfer. This encourages the User to remain active as long as possible. It easily transfers Users from bed to chair, chair to commode/bathroom etc. in a fraction of the time required for a conventional sling lifter.

- » Low castors for better access under beds
- » Large knee pads for safe and secure transfers
- » Dewert Control System with emergency Stop button
- » Ergonomic push handle
- High weight capacity of 175 kg (27.5st)

Dewert Control System with off board charging system







#### Stand Up Sling\*

The Stand up sling is designed to assist with the standing process. Fully padded polyester sling with non slip mesh and adjustable body belt for added security. Only suitable for those who have a degree of weight bearing capability

\*Slings fit all Parklite and Novaltis Hoists

CODE	DESCRIPTION
3901BE	3901BE Small Stand Up Sling
3902BE	3902BE Medium Stand Sling
3903BE	3903BE Large Stand Up Sling

### Hermes Hoist

The Hermes 250 is a very robust bariatric hoist. This bariatric hoist is fitted with larger 10cm castors, making manoeuvrability of the hoist easier. This hoist has a Linak Control System.

The battery is charged on board and the mains lead is neatly stored on a mains lead retainer on the mast of the hoist. The four point spreader bar offers greater comfort for the Client when being hoisted in a sling. The spreader bar is also angled to offer greater head clearance. There are metal clips fitted to each spreader bar point, ensuring no movement of the sling loops when fitted.

#### Portable Bed Hoist



				5	
<b>CODE</b> 18011BE	LIFTING RANGE (BOOM POSITION LOW-HIGH) 54 - 180cm (21 - 71″)		CHASSIS WIDTH (INTERNAL) 60 - 110.5cm (23 - 43")	CHASSIS HEIGHT 11.4cm (4″)	WEIGHT CAPACIT 250kg (39st)

CODE	DESCRIPTION	HEIGHT (FLOOR TO TOP OF BOOM)	OVERALL WIDTH	WIDTH WITH EXTRA
HST13	Overhead Hoist	204.8cm	202cm	283.5cm

197.5cm

50 Kg

62 Kg



#### Slings

When choosing a sling, consider the following factors:- the amount of support required, the tasks which need to be undertaken, the comfort of the person being lifted, the ability of the person being lifted, the ability of the Carer and the hoist used.

- All Drive DeVilbiss slings have colour coded edges so that you can easily determine its size:- Black = Extra Large, Green = Large, Yellow = Medium, Red = Small
- Each sling has colour coded positional loops to offer different seating positions
- All slings come with padded legs as standard for optimum comfort. The padded legs are fitted with nylon inserts to stop the legs from bunching
- » A modesty strap is fitted to the left hand leg of the Fast Fit Deluxe Sling
- Nylon inserts are fitted to the head area of the Comfort Sling for additional support
- The rear of the slings have straps fitted to assist the Carer to position and apply the sling to the client in safety
- Slings manufactured from spacer fabric are also available. Spacer fabric is a 3 dimensional material, which is more suitable for sitting on for longer periods of time
- All slings have woven labels for greater longevity and have a QR code, when scanned, takes you through to a video of how to assemble the slings
- Maximum user weight limit of slings:- Nett (200Kg), Spacer (300Kg) Polyester (400Kg)

SLING	EXTRA SMALL	SMALL	MEDIUM	LARGE
Fast Fit Sling	3201XBE	3201BE	3202BE	3203BE
Deluxe Hammock Sling	4501XBE	4501BE	4502BE	4503BE
Fast Fit Deluxe Sling - Polyester	3401XBE	3401BE	3402BE	3403BE
Fast Fit Deluxe Sling - Net	3404XBE	3404BE	3405BE	3406BE
Comfort Sling	3301XBE	3301BE	3302BE	3303BE
Dress Sling	3101XBE	3101BE	3102BE	3103BE
Split Leg in Chair Hammock Sling	4001XBE	4001BE	4002BE	4003BE
In Chair Hammock Comfort Sling	4101XBE	4101BE	4102BE	4103BE
High Back Toilet Sling	4201XBE	4201BE	4202BE	4203BE



4003XBE
4103XBE
4203XBE

The Fast Fit Sling in polyester is a universal all purpose sling, fast and easy to apply with full back and leg support. The padded legs are fitted with nylon inserts to stop the legs from bunching. This sling can be fitted to cross the Users legs to protect their modesty. Secured by loop fastening.





The Deluxe Hammock Sling is a hybrid sling designed to offer the full support of a hammock sling with the function of a universal. Can be applied and removed while seated in a chair, whilst the full padded support of the upper thigh gives maximum comfort, lifting from both inner and outer thigh to reduce hip rotation and abduction. This sling is suitable for the majority of users including amputees and those prone to excessive movements. Secured by loop fastening.

> The Fast Fit Deluxe Sling is designed to give full body support, offering additional support to sacral and thigh area. It is available in both polyester and net fabric. A modesty strap is fitted to the left hand leg. Secured by loop fastening.





The Split Leg In Chair Hammock Sling comes with hip straps and pockets to tuck lifting straps into and is supplied in black spacer fabric. This sling has a split leg design allowing you to fold away the leg supports preventing build up of heat and adding to user comfort. The spacer fabric offers additional benefits which are the fabrics 3 dimensional and used for the replacement of foam padding. The fabric is breathable, very durable and provides excellent pressure distribution. Secured by loop fastening.

The In Chair Hammock Comfort Sling was designed to give a more comfortable transfer for Users who experience discomfort when hoisted. A seamless bottom edge reduces pressure caused by taught binding and is complimented with a thinner padded pommel. This sling has the benefit of spacer fabric, reducing heat build-up and through 4 way stretch as well as being seamless, helping with pressure distribution. Suitable for those who are unable to remove slings after transfer.



The Comfort Sling in polyester has been designed to give extra support and comfort. It is a general purpose sling with an integral head support. Nylon inserts are fitted to head area for additional support. Secured by loop fastening.

The High Back Toilet Sling provides additional support, enabling those users who can support themselves sufficiently enough to resist the gravitational effect of a large aperture, to safely toilet directly from a seated position. With a wide body belt, high and rigid back and the use of extensor straps anchored to give maximum support along with a unique low cut around the hips, the problem of under arm pressure is eliminated.



The Dress Sling in polyester is an easy fit sling, specifically designed for easy access for toileting/personal hygiene. This sling allows access to adjust clothing and is fitted to head area for additional support.





#### **Acute Riser Recliner**

Ideal for individuals who have difficulty standing up and sitting down due to restricted mobility e.g. those with MS, Arthritis, Rheumatism or circulation problems. The Acute Riser Recliner Armchair has been specifically designed for the Acute Sector and can be easily manoeuvred from one side of the room to the other.

- » One motor activates a distinct lift, rest and recline action
- » It will rise and tilt users to their feet
- » Quiet motors provide a smooth and discreet operation
- » Ergonomically contoured seat and back for additional comfort
- » Fully padded chaise for complete comfort
- » Comes with a simple to operate handset
- » 2 year warranty on wooden frame and electrical components
- » Anti-microbial and crib 5 fire retardant PVC fabric

**Colour Options** 

Black PVC

Press down on pedal to lower wheels to manoeuvre the chair with ease

### **Chicago Rise & Recliner**

The Anti-Microbial Chicago has been tested to BS 5852 ignition source 5 (Crib 5) making it ideal for nursing and care home environments. The easy to clean PVC has passed the Anti-Microbial standard ASTM G-21-96.

- This Chicago model has been Crib 5 tested, so it is suitable for Nursing, Care and Residential Homes
- Anti-Microbial PVC fabric helps minimise build up of harmful bacteria such as MRSA and is easy to wipe clean for improved infection and hygiene control
- » Stylish and modern upholstery to suit any decor
- Ergonomically contoured seat and back and padded upholstery for added comfort and support
- » One motor activates a distinct lift, rest and recline action
- » It will rise and tilt Users to their feet
- Excellent snooze position with backrest reclining to 64° for supreme comfort and relaxation
- » Quiet motors provide a smooth and discreet operation
- » Quality assured, meeting all relevant fire retardancy standards
- Safety key is built into hand controller to prevent accidental movement and entrapment
- » Stitched handset holder to prevent loss or damage to the handset
- Ideal for individuals who have difficulty standing up and sitting down due to restricted mobility e.g. those with MS, Arthritis, Rheumatism, circulation problems
- » Dual handset model available as a special order (CLR15)

#### Colour Options





CODE	DESCRIPTION	WIDTH	HEIGHT	DEPTH	SEAT WIDTH	SEAT HEIGHT	SEAT DEPTH	WEIGHT CAPACITY
CLR11	Chicago AM PVC Fabric	83cm (33″)	108cm (43")	93cm (37″)	50cm (20″)	48cm (19″)	52cm (20″)	150kg (23st)

CODE	DESCRIPTION	WIDTH	HEIGHT	DEPTH	SEAT WIDTH	SEAT HEIGHT	SEAT DEPTH	WEIGHT CAPACITY
CLR27BLK	Acute Riser Recliner	82cm (32")	105cm (41")	94cm (37")	53cm (20")	48cm (19″)	55cm (21″)	130kg (20st)

Memory Foam Seat

#### **Integra Shell Seat**

The Integra Shell Seat tilts forward easily to help the User stand up from the chair. The chair is useful for Users who are seated for long periods at a time and at risk of developing pressure ulcers, therefore requiring additional support, comfort and positioning assistance.

- » Tilt in space chair gas strut operated
- » The chair will rise and tilt Users to their feet
- » Easy inclination due to the ergonomically designed push handle with hand lever
- » Angle adjustable and removable legrest which automatically moves under the seat when the chair is tilted forwards, enabling the User to get out of the chair with ease
- » Velcro attached ergonomically designed headrest, which can be easily adjusted
- » Lap belt is fitted as standard for additional safety
- » Handy storage pocket on the backrest
- » Easy to manoeuvre high quality brake wheels (12.5cm) which integrally rotate 360 degrees
- » Directional wheel enables the chair to be transferred in a straight line
- » Good hoist access for easy moving and handling
- » Brown leather style upholstery on the outside, together with 2 way stretch, breathable Dartex material on the inside
- » High risk cushion supplied as standard
- » Optional footrest & table available
- » Integra chair available for the Acute Sector c/w waterproof seat cushion (code SHELLINTACBE)

			5.0				
CODE	SEAT WIDTH	SEAT DEPTH	SEAT HEIGHT	OVERALL WIDTH	OVERALL DEPTH	ARMREST HEIGHT	WEIGHT CAPACITY
SHELLINT16BE	39cm (15.5″)	41cm (16")	45cm (18″)	67 - 72cm (26 - 28″)	89cm (35″)	23cm (9″)	120kg (19st)
SHELLINT18BE	45cm (18″)	44cm (17")	49cm (19″)	67 - 72cm (26 - 28″)	89cm (35″)	23cm (9″)	120kg (19st)
SHELLINT20BE	51cm (20″)	47cm (18.5″)	53cm (21″)	67 - 72cm (26 - 28″)	89cm (35″)	23cm (9″)	120kg (19st)

Ergonomic wing

headrest

inclination thanks shape and the

CODE	DESCRIPTION
SHELLINT16	Integra Chair - 16"
SHELLINT18	Integra Chair - 18"
SHELLINT20	Integra Chair - 20"
SHELLINT18AC	Integra Chair - 18″ - Acute

### Mountway Solo Toilet Lift

The Solo Toilet Lift replicates natural standing movement at the touch of a button, assisting Users to help themselves get up into a standing position.

- Battery powered device that enables Users with restricted mobility to get on and off toilet at the touch of a button
- Ideal for individuals who have difficulty performing sit to stand activities due to muscle or joint weakness in their legs and/or upper body e.g. due to multiple sclerosis, osteo or rheumatoid arthritis, hip replacement, heart problems and obesity
- Helps reduce manual handling risks for Carers who provide toileting assistance
- » Powered lifting action replicates body's natural standing movement
- Seat can be stopped at any position to suit User comfort and support requirements
- Fits easily and compactly over existing toilet, blending discreetly into modern home
- » Smooth, wipe clean surfaces for ease of cleaning
- Acts as a normal toilet seat in lowered position for use by other household members
- » No specialist installation simply remove existing toilet seat and place over toilet bowl
- Feet rotate clockwise and anti clockwise to adjust height for different Users/toilets
- Models also available with tilt up armrests (with/without integral operating controls), which give additional User support
- Supplied with two rechargeable batteries to ensure continuous operation
- » Battery level indicator monitors when recharging is required
- Set of extended feet is available as an accessory to raise seat a further 6cm (2.5") for accommodating tall Users or high toilets (Product code: TL634)

CODE	DESCRIPTION
MSOHSUK	Solo Toilet Lift without Arms
MS2LSUK	Solo Toilet Lift with 2 Arms (Left Switch)
MS2RSUK	Solo Toilet Lift with 2 Arms (Right Switch)

Models also available with tilt up armrests (with, without integral operating controls which give addition User support

CODE	WIDTH	DEPTH	HEIGHT	SEAT WIDTH	SEAT HEIGHT DOWN	SEAT HEIGHT UP	CLEARANCE HEIGHT	CLEARANCE WIDTH	WEIGHT CAPACITY
MSOHSUK	61cm (24")	54cm (21")	43cm – 49cm (17″ – 19″)	55cm (22")	43cm – 49cm (17″ – 19″)	56cm - 62cm (22" - 24")	36cm – 42cm (14″ – 17″)	39cm (15″)	160kg (25st)
MS2LSUK	71cm (28")	54cm (21")	66cm – 72cm (26″ – 28″)	55cm (22")	43cm – 49cm (17″ – 19″)	56cm - 62cm (22" - 24")	36cm – 42cm (14″ – 17″)	39cm (15″)	160kg (25st)
MS2RSUK	71cm (28″)	54cm (21")	66cm – 72cm (26″ – 28″)	55cm (22")	43cm – 49cm (17″ – 19″)	56cm – 62cm (22″ – 24″)	36cm – 42cm (14" – 17")	39cm (15″)	160kg (25st)




#### **AquaJoy Premier Plus**

There are a wide range of bathlifts and accessories available. For more please visit www.drivedevilbiss.co.uk

The Premier Plus Bathlift is a market leading, reclining, versatile bathlift that fits almost any style of bath including corner baths. It is suitable for both adults and children with a high weight capacity and a wide range of accessories.


- This reclining bathlift allows users to lower / raise themselves safely and smoothly in and out of the bath
- Reclining angle up to 40° makes full use of bath space to maximise legroom for transfers / bathing
- Separates easily into two lightweight sections for easy handling – backrest simply folds forward to detach
- High seat travel allows safe, level transfers in deeper baths and low travel offers excellent water immersion
- Sturdy, plastic construction has highest weight capacity of any reclining bathlift – ideal for heavier users up to 170kg (27st)
- Supplied with side flap protectors to allow smooth movement over bath handles
- Lifting bar on easy release suction feet assists to remove bathlift with ease

- Can be used with or without cover mats; these are easily removed for machine washing
- Lightweight, floating, fully waterproof hand control is comfortable to hold with soft touch, easy to operate buttons for users with limited dexterity or sight impairment
- Battery indicator light signals red when recharging is required
- Built-in safety feature prevents the bathlift from lowering unless there is sufficient charge to raise users to the top of the bath
- Wide range of accessories offer varying levels of support for both adults and children
- Suitable for thermal decontamination in automated washing machines (excluding battery)

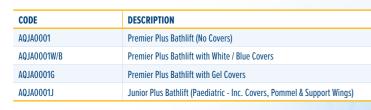

#### PREMIER PLUS TECHNICAL DATA

HEIGHT MINIMUM	9cm (3.5″)	BACKREST HEIGHT	60cm (24″)
HEIGHT MAXIMUM	46cm (18")	BACKREST WIDTH	34cm (13")
WIDTH (INC. SIDE FLAPS)	69cm (27")	BACKREST RECLINE	40°
DEPTH (WITH BACKREST RECLINED)	90cm (35″)	BASEPLATE WIDTH	39cm (15″)
DEPTH (WITH BACKREST UPRIGHT)	65cm (26")	BASEPLATE DEPTH	58cm (23″)
SEAT WIDTH	39cm (15″)	WEIGHT (WITHOUT BATTERY)	13kg (29lb)
SEAT DEPTH	49cm (19")	MAXIMUM WEIGHT CAPACITY	170kg (27st)

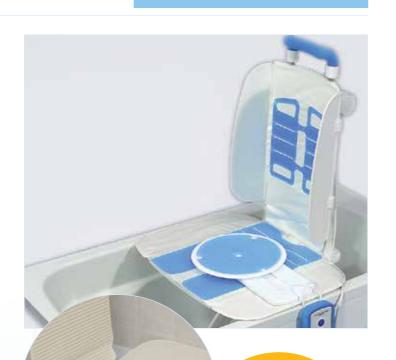

#### Swivel & Slide

sists users to transfer on and off bathlift. Features gel filled padding for comfort

#### **Mountway Neptune**

The Neptune Bathlift is one of the lightest battery powered fixed back bathlifts available, and is suitable for both corner baths and wider baths for maximum flexibility.

- Powered seat gently and safely lowers and raises Users smoothly in the bath at the touch of a button to maintain hygiene and soothe aching joints in warm water
- Compact design allows the bathlift to fit into almost any style of bath for optimal flexibility
- Slim line structure sits at the very back of the bath giving impressive legroom and allowing space for transfers / bathing later of the structure for the structure for the structure of the structur
- Large 13cm (5") suction feet grip securely to bath surface for added safety
   Contoured side flaps avoid most bath handles for unobstructed use
- Sits at only 7cm (3") from the bottom of the bath allowing improved water immersion
- Optional extra wide side flaps provide a safe transfer platform and maximise legroom in corner baths
- One of the lightest fixed back bathlifts available at 9.5kg (21lb)
- » Compact design allows for space saving storage
- » Separates easily into 2 lightweight sections for quick fitting and removal
- Smooth and accessible surfaces enable quick and simple cleaning of the lift and the bath
- Antimicrobial solution in suction feet, base frame and handset assists with infection and hygiene control
- Waterproof, slim line hand control with large textured rocker switches can be operated by people with limited dexterity or visual impairment
- » Battery level indicator light monitors when recharging is required and Stop Descent<sup>™</sup> feature prevents bathlift from lowering without sufficient battery power to rise
- Seat is made from sturdy glass filled polypropylene (PP) with a robust, rust proof aluminium frame
- Suitable for thermal decontamination in washing machines (excluding battery)
- Warranty 5 years on frame/seat, 2 years on consumables and 1 year on battery and accessories

NEPTUNE TECHNICAL DATA			
HEIGHT MINIMUM	7cm (3″)	BACKREST HEIGHT	66cm (26″)
HEIGHT MAXIMUM	43.5cm (17")	BACKREST WIDTH	34cm (13")
WIDTH (INC. SIDE FLAPS)	73cm (29″)	BASEPLATE WIDTH	33.5cm (13″)
DEPTH (UPRIGHT)	58cm (23″)	BASEPLATE DEPTH	39cm (15″)
SEAT WIDTH	35cm (14")	WEIGHT (WITHOUT BATTERY)	9.5kg (21lb)
SEAT DEPTH	41cm (16″)	MAXIMUM WEIGHT CAPACITY	160kg (25st)

There are a wide range of bathlifts and accessories available. For more please visit www.drivedevilbiss.co.uk



CODE	DESCRIPTION
NEPA	Neptune Bathlift

## **Bathroom Accessories**

### Vitaturn K - Turning Aid

Allows Users with limited mobility to turn on seat for easier transfers into bath. Smooth rotating operation for added comfort. 2cm (1") height on seat still maintains good water immersion.



2322

MILLING

CODE	DESCRIPTION
412100100	Vitaturn K - Turning Aid

#### Vitaturn L - Turning Aid & Transfer Aid

Securely attaches to the sliding surface offering more safety and relief for the Carer. Assists Users with restricted mobility to slide onto seat from outside bath and turn to lift their legs into bath.



CODE	DESCRIPTION
101385	Vitaturn L - Turning and Transfer Aid





#### Leaflet Code: LL222

For more information contact your local dealer:



#### www.drivedevilbiss.co.uk

All sizes and weights stated are nominal. All errors and omissions excepted. Due to the limitations of photography and the printing process colours shown may not be 100% accurate.